




# Harmar

AMERICA'S LIFT LEADER

## IVCLIFT10 & IVCLIFT14 Invacare® Exclusive Power Wheelchair Lifts



Lift folds automatically  
when not in use

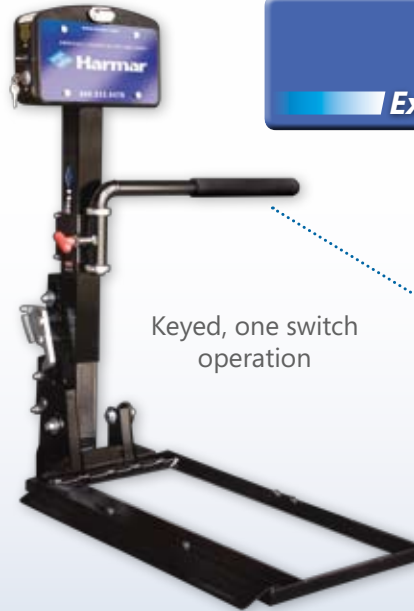

INVACARE  
Exclusive Product

Keyed, one switch  
operation



Adjustable padded hold  
down arm provides a secure,  
custom fit for every chair  
without the use of straps or  
tie-downs

Lightweight, unobtrusive platform designed specifically for use with  
Invacare® Pronto®\*\* and TDX®\*\*\* series power wheelchairs

### FEATURES

- Designed exclusively for Invacare
- Padded hold-down arm conveniently secures power wheelchair
- No power wheelchair modifications required
- Lift folds automatically when unoccupied
- No straps or tie-downs required
- Manual crank backup & license plate mount included
- Swing-away option available

### SPECIFICATIONS

Lifting Capacity:	IVCLIFT10 - 300 lbs IVCLIFT14 - 360 lbs
Platform Size:	7.5" x 27"
Installed Weight:	89 lbs
Hitch Height (range):	5"-24"*
Hitch Class:	IVCLIFT10 - Class II or III IVCLIFT14 - Class III
Warranty:	3-Year Transferable

\*With appropriate hitch adapter installed

### DESCRIPTION

**Compact, durable and American-made, the IVCLIFT10 and IVCLIFT14 outside vehicle platform lifts safely and securely transport Invacare® Pronto®\*\* and TDX®\*\*\* series power wheelchairs. With a tidy and space efficient rear profile, this lift has a highly desirable and understated presence on the vehicle exterior.**

The IVCLIFT10 and IVCLIFT14 are easy to operate and are built to last for years of trouble-free use. Designed to meet the needs of consumers who need a quick, automatic securing option without making modifications to their power wheelchair, the stylish and compact detailing of this lift make it a must have for Invacare TDX and Pronto series power wheelchair owners.

\*\* IVCLIFT10: Invacare Pronto M41, M51 & M61 power wheelchairs

\*\*\* IVCLIFT14: Invacare TDX-SI, TDX-SP & TDX-SR power wheelchairs. TDX-SR models with non-powered seating ONLY.

© 2011 Harmar. Invacare, The Blue Medallion Design, TDX and Pronto are registered trademarks of Invacare Corporation. M41, M51 and M61 are trademarks of Invacare Corporations.  
Form No. 11-030

800-833-0478 • www.harmar.com



# Harmar

AMERICA'S LIFT LEADER

## IVCLIFT10 & IVCLIFT14 Invacare® Exclusive Power Wheelchair Lifts



### IVCLIFT10 Power Wheelchair Lift

For transport of Invacare® Pronto® M41™, M51™ and M61™ power wheelchairs



### IVCLIFT14 Power Wheelchair Lift


For transport of Invacare® TDX®-SI, TDX®-SP and TDX®-SR power wheelchairs

### Easy as 1, 2, 3!



1

Drive chair onto platform

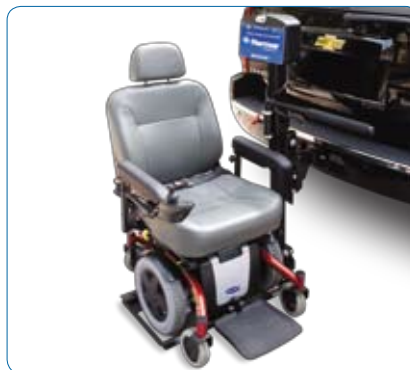

2

Raise armrest closest to vehicle and rotate swing arm out over top of power wheelchair



3

Raise lift and go!

800-833-0478 • [www.harmar.com](http://www.harmar.com)