

# TRIFOLD® Advantage Series® RAMP



The TRIFOLD® Advantage Series® portable wheelchair ramp, with its unique 3-fold design, offers the length required for wheelchairs and scooters to easily access steps, vehicles, and raised landings. The ramp is designed to be used and carried as one unit, or if desired, can be separated into two individual sections.

Innovative hinge design features interlocking brackets to provide smoother operation and reduced pinch points. Ramp features an applied skid-resistant surface with reflective safety stripes, as well as self-adjusting bottom transition plates for easy conversion from ramp to ground. Manufactured of welded aircraft-grade aluminum.

Patented.




Bottom transition plates adjust to uneven terrain



Ergonomically designed handles offer convenient carrying and a comfortable grip




Separates for transport and storage

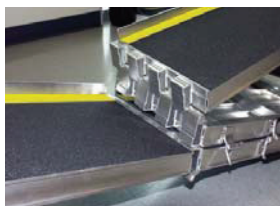

Hinge

## OPTIONAL TOP LIP EXTENSION (TLE)

Extend the usefulness of your TRIFOLD® Advantage Series® portable wheelchair ramp with the Top Lip Extension (TLE). This optional plate increases the standard 3" lip to become 9", allowing the ramp to clear the rear bumper of SUVs and vans by simply attaching to the top lip. Shipped as a pair.

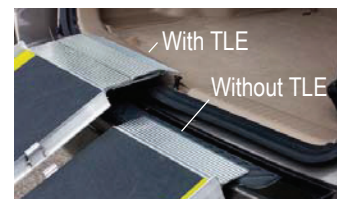

SPECIFICATIONS	TRIFOLD AS5	TRIFOLD AS6	TRIFOLD AS7	TRIFOLD AS8	TRIFOLD AS10	TLE
Made of aluminum	•	•	•	•	•	•
Non-skid surface	•	•	•	•	•	•
Usable size	5' x 29½"	6' x 29½"	7' x 29½"	8' x 29½"	10' x 29½"	8½" x 137/8" x 23/8"
Folded size (L x W x H)	30¾" x 15" x 11½"	36¾" x 15" x 11½"	42¾" x 15" x 11½"	48¾" x 15" x 11½"	60¾" x 15" x 11½"	n/a
Built-in carrying handle	•	•	•	•	•	n/a
Weight (lbs.) *	16 ea. (2 sections)	19 ea. (2 sections)	22 ea. (2 sections)	25 ea. (2 sections)	30 ea. (2 sections)	5 per set
Weight capacity (lbs.)	800	800	800	800	800	**500
For wheelchairs & scooters	•	•	•	•	•	•
Usage guidelines	Curbs, 1-2 steps, some minivans	1-3 steps, minivans, some vans & SUVs	1-3 steps, minivans, some vans & SUVs	2-3 steps, minivans, some vans & SUVs	2-4 steps, minivans, some vans & SUVs	n/a

\* Weights are approximate and may vary. Weights listed are per ramp section (each section can be carried, assembled, and stored separately). Consult your equipment's owner guide for proper degree of incline. Never exceed its recommendations.

\*\* Ramp weight capacity reduced when used with the TLE.



Bariatric weight capacity = 800 lbs. or greater



UPS Shippable